

SOMNI
SCIENTIFIC



SOMNI REBREATHING ANAESTHESIA SYSTEM



QUALITY:

- H-frame provides exceptional support and stability.
- 3" wheels for ease of transport.
- High quality, durable, limited maintenance design.
- Silicone tubing and gaskets throughout.
- Designed to switch from a rebreathing circuit to a non-rebreathing circuit with a simple toggle switch.

▶ COMPATIBLE ACCESSORY OPTIONS

Dual Vaporiser
Bain Block Manifold
Dual Gas supply
Dual E-cylinder
O2 Concentrator mount
5 year warranty




SAFETY:

- Pressure Manometer center mount for immediate visibility.
- Popoff valve releases excess pressure at 20cm H₂O while maintaining passive volume in the rebreathing bag.
- The oxygen toggle switch works as a built-in safety system, that keeps the pressure at 20PSI. This is designed with patient and user safety in mind.




CONVENIENT:

- Lightweight, small footprint, easy assembly.
- Easily Configured shelves for additional workspace or accessories.
- Powder coated and anodized frame for easy cleaning.

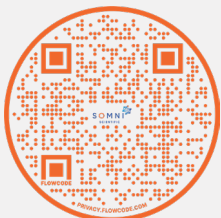

(T) 0800 0129101 (D) 01872 248890 (M) 07798 969805

enquiries@somniscientific.com  www.somniscientific.co.uk