

# SOMNI

SCIENTIFIC



## INDUCTION CHAMBERS ACTIVE VACUUM AND PASSIVE STYLES

### Induction Chambers

#### Active Vacuum and Passive Styles

The SOMNI line of induction chambers are designed with slide top lids that allow for one handed operation. Unlike a “hinge lid” the slide top keeps gases from being “drawn” into the workplace when the chamber is opened.


### ► SPECIFICATIONS

**ACTIVE LARGE:** 2.31 kg / 29.21 cm L 15.24 W 15.24 H

**ACTIVE MEDIUM:** 1.4 kg / 23.5 cm L 12.7 cm W 12.7 cm H

**PASSIVE MEDIUM:** 1.4 kg / 23.5 cm L 12.7 cm W 12.7 cm H

**PASSIVE LARGE:** 2.31 kg / 29.21 cm L 15.24 W 15.24 H

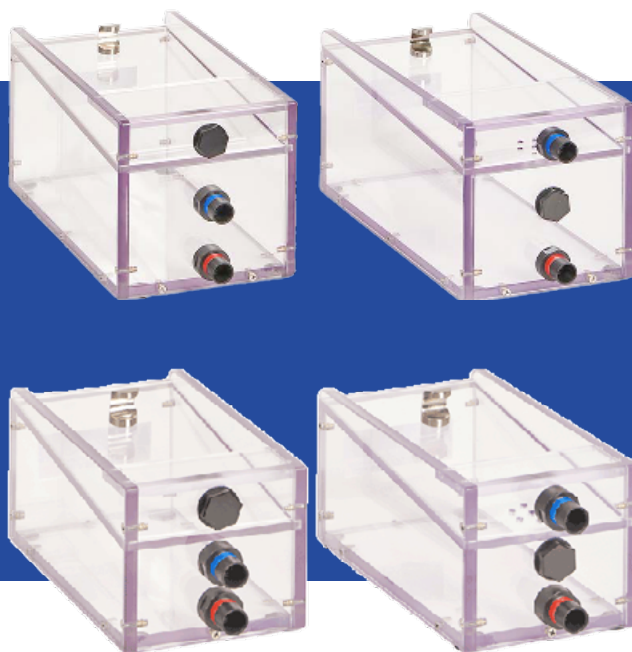

### QUALITY:



- Constructed of durable, shatter resistant polycarbonate.
- Can be safely cleaned with common clinical disinfectants.

### SAFETY:

- Prevents unnecessary exposure to personnel and environment (NIOSH 2ppm) when used with active waste gas removal system EPS-3 or W.A.G. filter (active chamber only).

### CONVENIENT:

- One-handed operation
- Color coded inlet and outlet ports.
  - Compatible with 19mm tubing & 15mm adapters
- MRI compatible

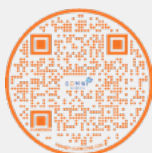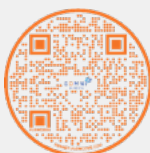

☎ (T) 0800 0129101 (D) 01872 248890 (M) 07798 969805

✉ enquiries@somniscientific.com 🌐 www.somniscientific.co.uk